

# New to range

## Orthopaedic GUIDE WIRES Single Use

The graduated spinal guide wires supplied by Judd Medical are used by various spinal Consultant Surgeons in complex procedures such as Multiple-Level Fusion and Correction of Scoliosis.


During surgery, the surgeon will locate the patient's pedicle, and once located, the Judd Medical spinal guide wires are put in place as a guide for the implantation of pedicle screws with the help of x-ray. The graduations on the wire provide the surgeon with information on the depth of the pedicle, making it easier for the placement of screws as they can calculate it based on how deep the wires are placed.

### Packaging

Sterile packaged for single use, saving on re-sterilisation costs

### Material

High quality Medical Grade 316L Stainless Steel

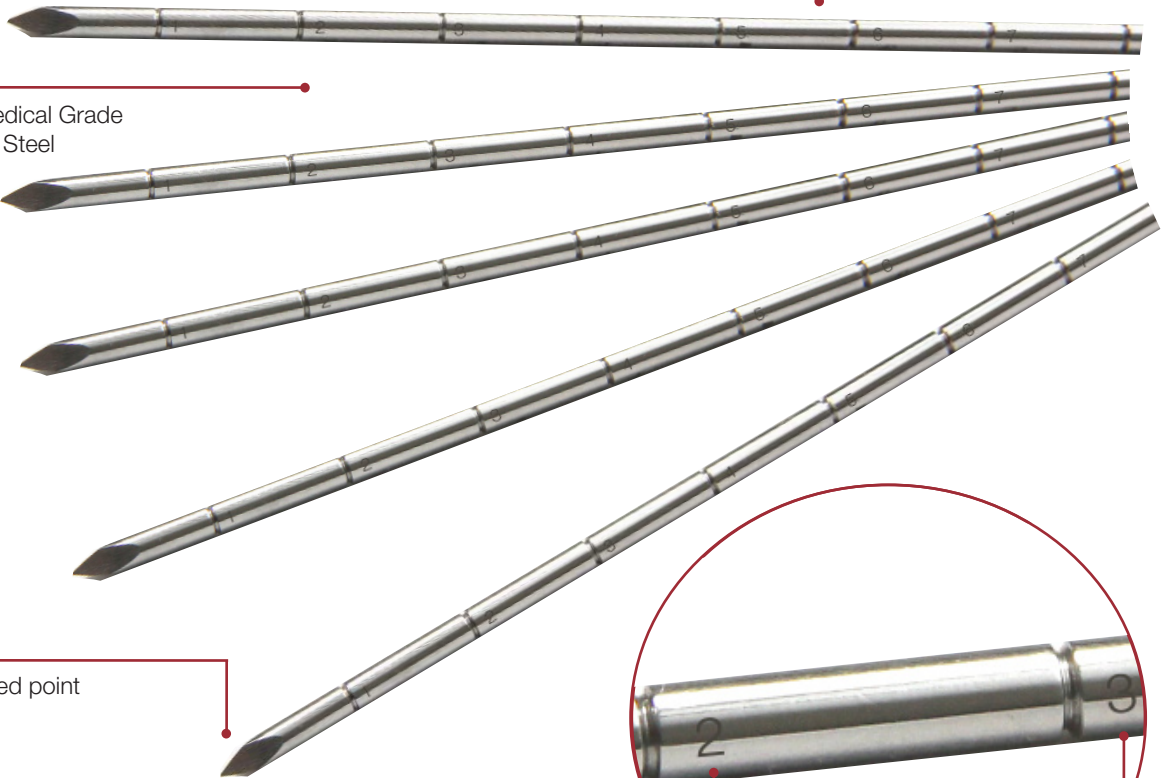


### Tip

Three-faceted point

### Accurately gauged

Indentations in increments of 1cm

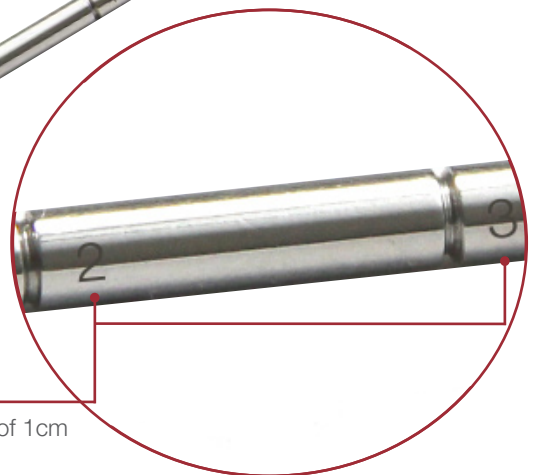

## Sterile Orthopaedic Guide Wires

Product Code	Description	Length	Single Use	Pk Size
<b>SG/2.0/18S</b>	Spinal Guide Wire 2.0mm	18cm	•	20