

SHUKLA™ Medical

Innovators of
Xtract-All®

Maxi Ratcheting Screwdriver System

Intended Use:

- Insertion and Extraction of Bone Screws

Includes:

- Solid Drivers
 - 8 Hex Drivers
 - 5 Torx Drivers
 - Cruciform Driver
 - Flat Tip Driver
 - Phillips Driver
- Flexible Drivers
 - 8 Hex Drivers
 - 5 Torx Drivers
 - Cruciform Driver
 - Flat Tip Driver
 - Phillips Driver



**Have the instrument you need
when you need it.**

Contact us at: **888-4shukla (888-474-8552)**
or **customer.service@ShuklaMedical.com**
See our website at: **www.ShuklaMedical.com**



BR-0001 Rev A

Maxi Ratcheting Screwdriver System

SHUKLATM
Medical

MKS1029		MAXI RATCHETING SCREWDRIVER SYSTEM		
Tip No.	Part Number	Description		QTY
1	MDS1123-M01	2.0mm Hex Driver (SOLID)		1
2	MDS1124-M01	2.5mm Hex Driver (SOLID)		1
3	MDS1125-M01	3.5mm Hex Driver (SOLID)		1
4	MDS1126-M01	4.0mm Hex Driver (SOLID)		1
5	MDS1127-M01	5/32" Hex Driver (SOLID)		1
6	MDS1128-M01	3/16" Hex Driver (SOLID)		1
7	MDS1129-M01	4.5mm Hex Driver (SOLID)		1
8	MDS1130-M01	5.0mm Hex Driver (SOLID)		1
9	MDS1131-M01	T-8 Torx Driver (SOLID)		1
10	MDS1132-M01	T-15 Torx Driver (SOLID)		1
11	MDS1133-M01	T-20 Torx Driver (SOLID)		1
12	MDS1134-M01	T-25 Torx Driver (SOLID)		1
13	MDS1135-M01	T-30 Torx Driver (SOLID)		1
14	MDS1136-M01	Cruciform Driver (SOLID)		1
15	MDS1137-M01	Flat Tip Driver (SOLID)		1
16	MDS1138-M01	Phillips Driver (SOLID)		1
1	MDS1139-M01	2.0mm Hex Driver (FLEXIBLE)		1
2	MDS1140-M01	2.5mm Hex Driver (FLEXIBLE)		1
3	MDS1141-M01	3.5mm Hex Driver (FLEXIBLE)		1
4	MDS1142-M01	4.0mm Hex Driver (FLEXIBLE)		1
5	MDS1143-M01	5/32" Hex Driver (FLEXIBLE)		1
6	MDS1144-M01	3/16" Hex Driver (FLEXIBLE)		1
7	MDS1145-M01	4.5mm Hex Driver (FLEXIBLE)		1
8	MDS1146-M01	5.0mm Hex Driver (FLEXIBLE)		1
9	MDS1147-M01	T-8 Torx Driver (FLEXIBLE)		1
10	MDS1148-M01	T-15 Torx Driver (FLEXIBLE)		1
11	MDS1149-M01	T-20 Torx Driver (FLEXIBLE)		1
12	MDS1150-M01	T-25 Torx Driver (FLEXIBLE)		1
13	MDS1151-M01	T-30 Torx Driver (FLEXIBLE)		1
14	MDS1152-M01	Cruciform Driver (FLEXIBLE)		1
15	MDS1153-M01	Flat Tip Driver (FLEXIBLE)		1
16	MDS1154-M01	Phillips Driver (FLEXIBLE)		1
—	MHRS0311	Cannulated Blue Silicone Ratcheting Screwdriver		1
—	CA174-M01	Case		1

Phillips® is a registered trademark of the Phillips Screw Company. Torx® is a registered trademark of Textron, Inc. Shukla Medical™ is not affiliated with the above mentioned companies and Shukla Medical™ makes no claim to copyrights or trademarks which are the property of the above mentioned companies. All rights of distribution, also by photocopy, reprint or storing and recovery in any kind of data processing are reserved and require our written approval. © Shukla Medical™ United States All rights reserved.



BR-0001 Rev A

Maxi Ratcheting Screwdriver System



Surgical Technique Guide

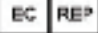




MKS 1029		Maxi Ratcheting Screwdriver System		
No.	Part Number	Description	Image	QTY
1	MDS1123-M01	2.0mm Hex Driver (SOLID)		1
2	MDS1124-M01	2.5mm Hex Driver (SOLID)		1
3	MDS1125-M01	3.5mm Hex Driver (SOLID)		1
4	MDS1126-M01	4.0mm Hex Driver (SOLID)		1
5	MDS1127-M01	5/32" Hex Driver (SOLID)		1
6	MDS1128-M01	3/16" Hex Driver (SOLID)		1
7	MDS1129-M01	4.5mm Hex Driver (SOLID)		1
8	MDS1130-M01	5.0mm Hex Driver (SOLID)		1
9	MDS1131-M01	T-8 Torx Driver (SOLID)		1
10	MDS1132-M01	T-15 Torx Driver (SOLID)		1
11	MDS1133-M01	T-20 Torx Driver (SOLID)		1
12	MDS1134-M01	T-25 Torx Driver (SOLID)		1
13	MDS1135-M01	T-30 Torx Driver (SOLID)		1
14	MDS1136-M01	Cruciform Driver (SOLID)		1
15	MDS1137-M01	Flat Tip Driver (SOLID)		1
16	MDS1138-M01	Phillips Driver (SOLID)		1
1	MDS1139-M01	2.0mm Hex Driver (FLEXIBLE)		1
2	MDS1140-M01	2.5mm Hex Driver (FLEXIBLE)		1
3	MDS1141-M01	3.5mm Hex Driver (FLEXIBLE)		1
4	MDS1142-M01	4.0mm Hex Driver (FLEXIBLE)		1
5	MDS1143-M01	5/32" Hex Driver (FLEXIBLE)		1
6	MSD1144-M01	3/16" Hex Driver (FLEXIBLE)		1
7	MDS1145-M01	4.5mm Hex Driver (FLEXIBLE)		1
8	MDS1146-M01	5.0mm Hex Driver (FLEXIBLE)		1
9	MDS1147-M01	T-8 Torx Driver (FLEXIBLE)		1
10	MDS1148-M01	T-15 Torx Driver (FLEXIBLE)		1
11	MDS1149-M01	T-20 Torx Driver (FLEXIBLE)		1
12	MDS1150-M01	T-25 Torx Driver (FLEXIBLE)		1
13	MDS1151-M01	T-30 Torx Driver (FLEXIBLE)		1
14	MDS1152-M01	Cruciform Driver (FLEXIBLE)		1
15	MDS1153-M01	Flat Tip Driver (FLEXIBLE)		1
16	MDS1154-M01	Phillips Driver (FLEXIBLE)		1
—	MHRS0311	Cannulated Blue Silicone Ratcheting Screwdriver		1
—	CA174-M01	Case		1

Contact us at: **888-4shukla (888-474-8552)**
 or **customer.service@ShuklaMedical.com**
 See our website at: **www.ShuklaMedical.com**

Maxi Ratcheting Screwdriver System

<p>Step 1 - Determine the drive configuration of the bone screw.</p>	<p>CRUCIFORM PHILLIPS TORX FLAT HEX</p> 
<p>Step 2 - Based on accessibility of the working site choose either the straight or flexible corresponding configuration screwdriver.</p>	
<p>Step 3 - Insert the 1/4" square end of the bone screwdriver into the blue silicon ratcheting handle. The driver shaft should click into place.</p>	
<p>Step 4 - Rotate the cap to the left for removing screws, the center position for locking the handle into place and to the right for installing screws.</p>	 <p>REVERSE LOCKED FORWARD</p>
<p>Step 5- The ratcheting screwdriver is now ready to remove or install bone screws.</p>	

				
<p>EMERGO GROUP Molenstraat 15 2513 BH, The Hague The Netherlands Tel: +31 (0)70 345 8570 Fax: +31 (0)70 346 7299</p>	<p>MKS1029</p>	<p>Shukla Medical 151 Old New Brunswick Road Piscataway, NJ 08854 U.S.A. 1(888)-4SHUKLA www.shuklamedical.com</p>	<p>CONSULT INSTRUCTIONS FOR USE</p>	<p>NON STERILE PRODUCT</p>